

# Ammonium 3

## Amonio

**Test 0-03**

REF 985 003



Rundküvettentest/ Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

0.04–2.30 mg/L  $\text{NH}_4\text{-N}$

0.05–3.00 mg/L  $\text{NH}_4^+$

0.05–3.00 mg/L  $\text{NH}_3$


690 nm

Method(e) / Método


**0031**

**0032**


**0033**

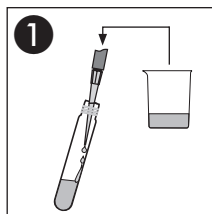


**4.0 mL** Probe  
Sample  
Echantillon  
Muestra

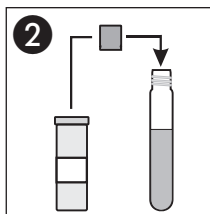


**1 x NANOFIX R2**

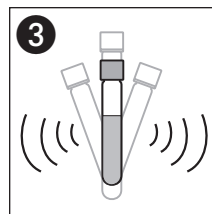


Schütteln  
Shake  
Agiter  
Agitar

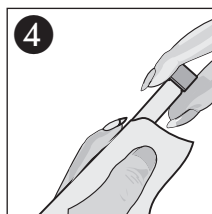

Säubern  
Clean  
Nettoyer  
Limpiar

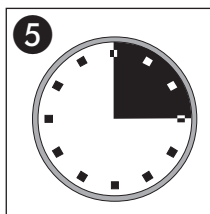

**15'00 min**

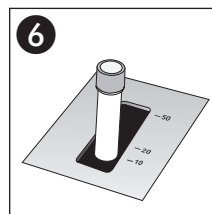

Messung  
Measurement  
Mesure  
Medición