

Ammonium 10

Amonio

Test 0-04

REF 985 004



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

0.2–8.0 mg/L $\text{NH}_4\text{-N}$

0.2–10.0 mg/L NH_4^+

0.2–10.0 mg/L NH_3


690 nm

Method(e) / Método


0041

0042


0043

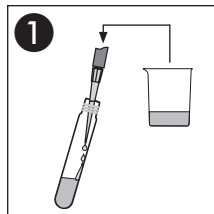


1
1.0 mL Probe
Sample
Echantillon
Muestra

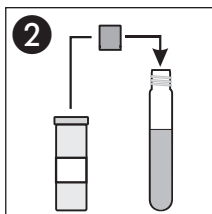


2
1 x **NANOFIX R2**

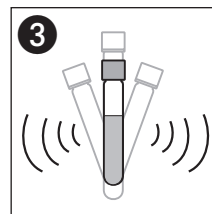


3
Schütteln
Shake
Agiter
Agitar

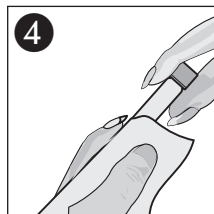

4
Säubern
Clean
Nettoyer
Limpiar

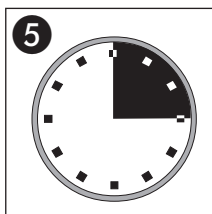

5
15'00 min

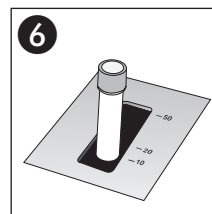

6
Messung
Measurement
Mesure
Medición