

# Ammonium 50

## Amonio

**Test 0-05**

REF 985 005



Rundküvettestest/ Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

1.0–40.0 mg/L  $\text{NH}_4\text{-N}$

1.0–50.0 mg/L  $\text{NH}_4^+$

1.0–50.0 mg/L  $\text{NH}_3$


690 nm

Method(e) / Método


**0051**

**0052**


**0053**

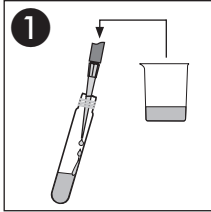


**200  $\mu\text{L}$**  Probe  
Sample  
Echantillon  
Muestra

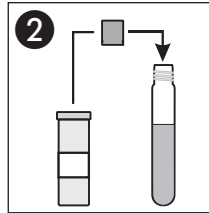


**1 x NANOFIX R2**

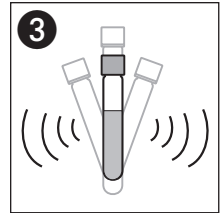


Schütteln  
Shake  
Agiter  
Agitar

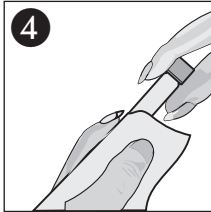

Säubern  
Clean  
Nettoyer  
Limpiar

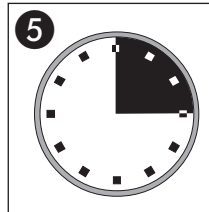

**15'00 min**

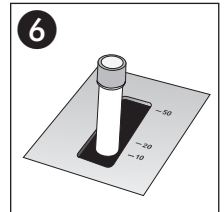

Messung  
Measurement  
Mesure  
Medición