

Ammonium 200

Amonio

Test 0-06

REF 985 006



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

30–160 mg/L $\text{NH}_4\text{-N}$

40–200 mg/L NH_4^+

40–200 mg/L NH_3


585 nm

Method(e) / Método


0061

0062


0063

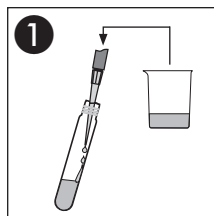


1
200 μL Probe
Sample
Echantillon
Muestra

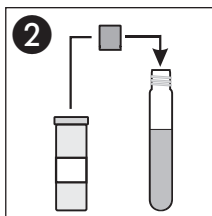


2
1 x **NANOFIX R2**

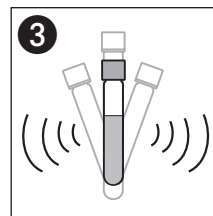


3
Schütteln
Shake
Agiter
Agitar

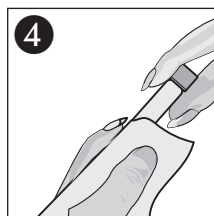

4
Säubern
Clean
Nettoyer
Limpiar

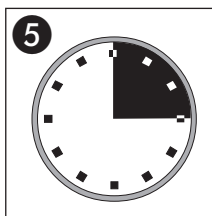

5
15'00 min

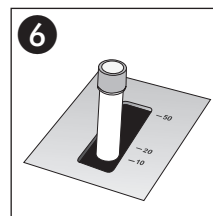

6
Messung
Measurement
Mesure
Medición