

# Ammonium 100

**Test 0-08**

REF 985 008



**Rundküvettentest/ Tube test**

**Test en cuves rondes / Prueba en cubeta cilíndrica**

4–80 mg/L  $\text{NH}_4\text{-N}$

5–100 mg/L  $\text{NH}_4^+$

5–100 mg/L  $\text{NH}_3$


**585 nm**

Method(e) / Método


**0081**

**0082**

**0083**

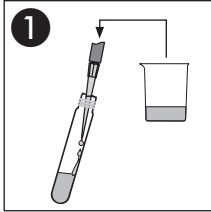


**500  $\mu\text{L}$**  Probe  
Sample  
Echantillon  
Muestra

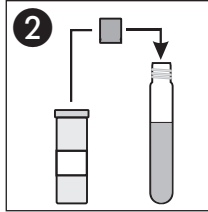


**1 x NANOFIX R2**




Schütteln  
Shake  
Agiter  
Agitar

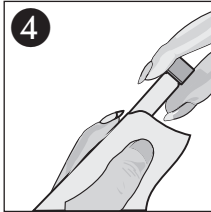

Säubern  
Clean  
Nettoyer  
Limpiar

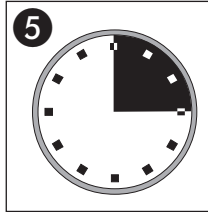

**15'00 min**

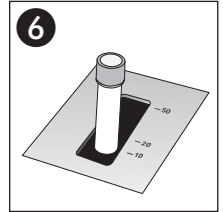

Messung  
Measurement  
Mesure  
Medición