

Blei 5

Lead/Plomb/Plomo

Test 0-09

REF 985 009



Rundküvettest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica


0.10–5.00 mg/L Pb

Messwert A / Sample value A / Valeur A / Valor de medición A:


520 nm

Method(e) / Método


0091

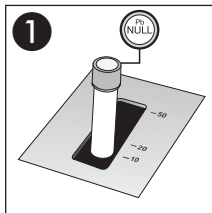


NULL messen
Measure blank
Lecture blanc
Lectura blanco

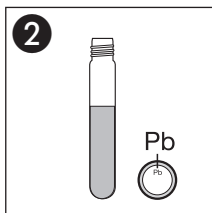


Blei-Rundküvette
Test tube **Lead**
Cuve ronde **Plomb**
Tubo de test **Plomo**

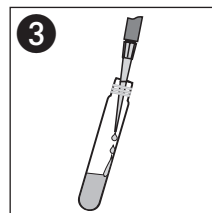


200 μ L R2

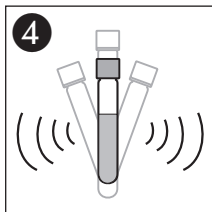


Schütteln
Shake
Agiter
Agitar

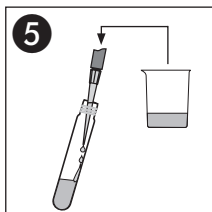


4.0 mL Probe
Sample
Echantillon
Muestra

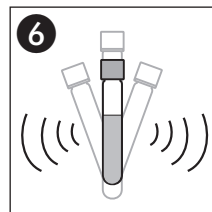


Schütteln
Shake
Agiter
Agitar

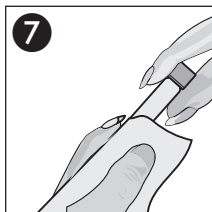

Säubern
Clean
Nettoyer
Limpiar

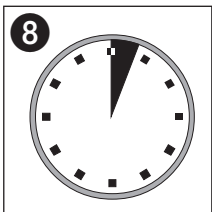

3'00 min

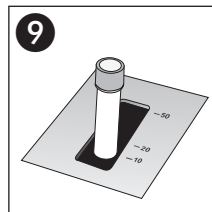

Messwert A
Value A
Valeur A
Valor A