

CSB 4000

COD/DCO/DQO

Test 0-11

REF 985 011



Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

400–4000 mg/L


620 nm

Method(e) / Método


0111

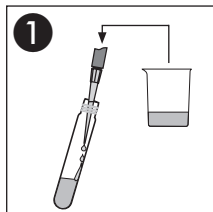


1.0 mL Probe
Sample
Echantillon
Muestra

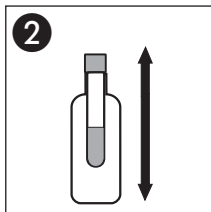


Sicherheitsgefäß
Safety bottle
Récipient de sécurité
Recipiente de seguridad

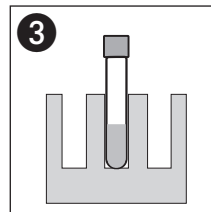


148 °C / 2 h

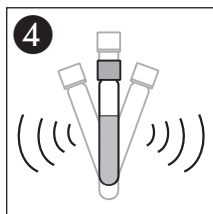


Schütteln
Shake
Agiter
Agitar

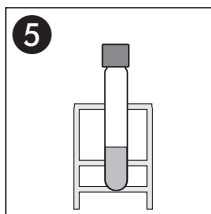

Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

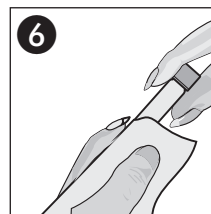

Säubern
Clean
Nettoyer
Limpiar

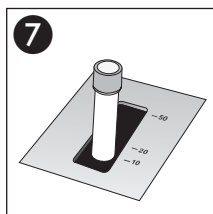

Messung
Measurement
Mesure
Medición