

CSB 60000

COD/DCO/DQO

Test 0-12

REF 985 012



Rundküvettest/ Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

5.0–60.0 g/L (5000–60000 mg/L)


620 nm

Method(e) / Método


0121

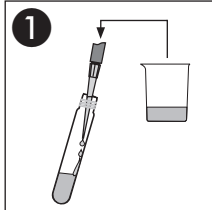


200 µL Probe
Sample
Echantillon
Muestra

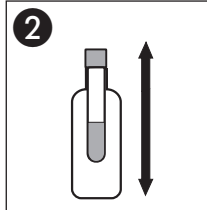


Sicherheitsgefäß
Safety bottle
Récipient de sécurité
Recipiente de seguridad

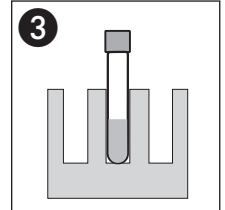


148 °C / 2 h

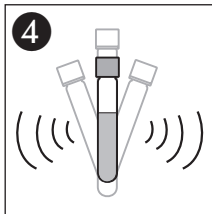


Schütteln
Shake
Agiter
Agitar

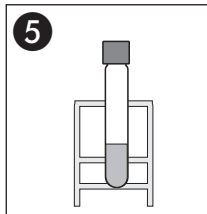

Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

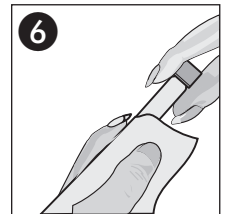

Säubern
Clean
Nettoyer
Limpiar

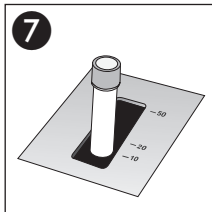

Messung
Measurement
Mesure
Medición