

Cadmium 2

Cadmio

Test 0-14

REF 985 014



Rundküvettentest / Tube test
Test en cuves rondes / Prueba en cubeta cilíndrica


520 nm

Method(e) / Método


0141

0.05–2.00 mg/L Cd²⁺


Nullwert / Blanc value / Zéro / Cero

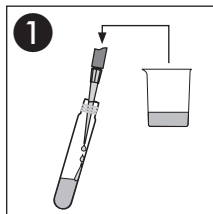


1
4.0 mL dest. Wasser
dist. water
d'eau distillée
agua dest.

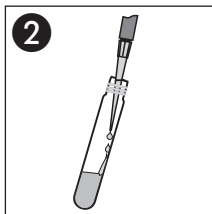


2
200 µL R2

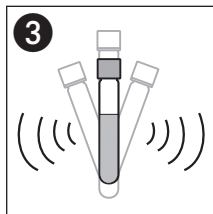


3
Schütteln
Shake
Agiter
Agitar

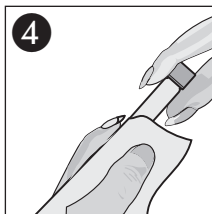

4
Säubern
Clean
Nettoyer
Limpiar

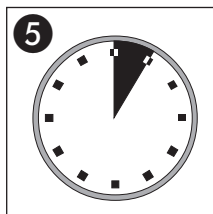


5
5'00 min

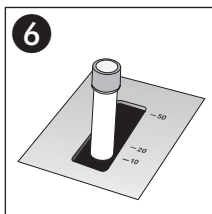


6
Messung
Measurement
Mesure
Medición


Messwert / Sample / Echantillon / Muestra

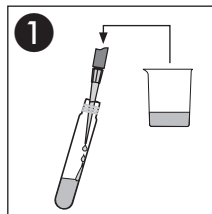


1
4.0 mL Probe
Sample
Echantillon
Muestra

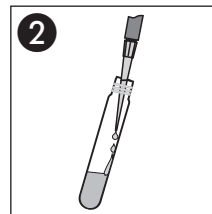


2
200 µL R2

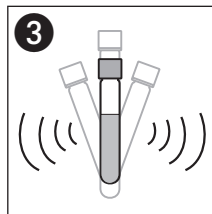


3
Schütteln
Shake
Agiter
Agitar

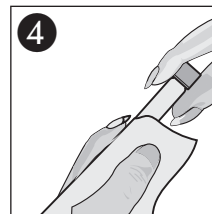

4
Säubern
Clean
Nettoyer
Limpiar

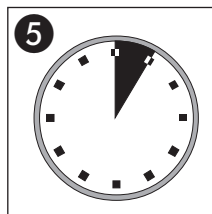

5
5'00 min

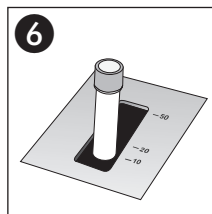

6
Messung
Measurement
Mesure
Medición