

Carbonathärte 15

Carbonate hardness / Dureté carbonatée

Alcalinidad de Carbonato

Test 0-15

REF 985 015



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

436/585 nm


Method(e) / Método

0151/0152


0153/0154

1.0–15.0 °d / 0.4–5.4 mmol/L H⁺


1.0–18.0 °e / 2.0–26.0 °f

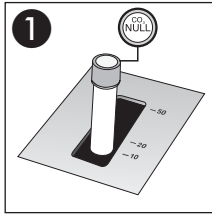


NULL messen
Measure blank
Lecture blanc
Lectura blanco

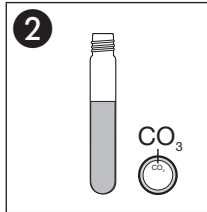


CO₃-Rundküvette
CO₃ test tube
Cuve ronde CO₃
Tubo de test de CO₃

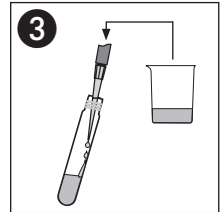


4.0 mL Probe
Sample
Echantillon
Muestra

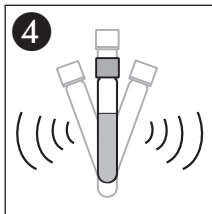


Schütteln
Shake
Agiter
Agitar

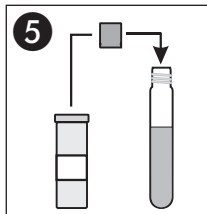


1 x NANOFIX R2

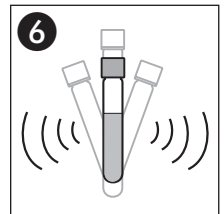


Kräftig schütteln
Shake well
Bien agiter
Agitar intensamente

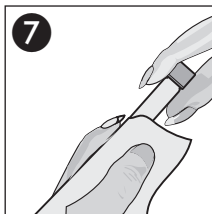

Säubern
Clean
Nettoyer
Limpiar

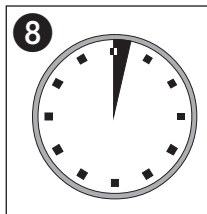

2'00 min

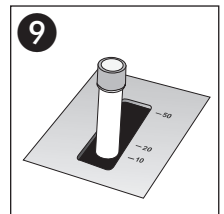

Messung
Measurement
Mesure
Medición