

Chlorid 50

Chloride/Chlorure/Cloruro

Test 0-21

REF 985 021



Rundküvettest/ Tube test


Test en cuves rondes/ Prueba en cubeta cilíndrica

0.5–50.0 mg/L Cl⁻


470 nm

Method(e) / Método


0211

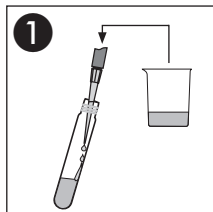


1
4.0 mL Probe
Sample
Echantillon
Muestra

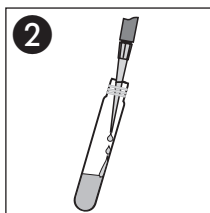


2
1.0 mL R2

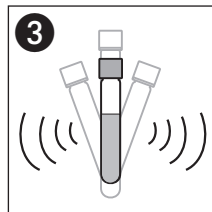


3
Schütteln
Shake
Agiter
Agitar

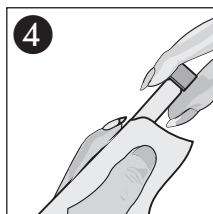

4
Säubern
Clean
Nettoyer
Limpiar

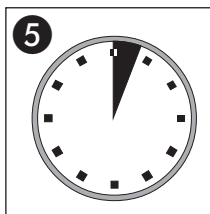

5
3'00 min

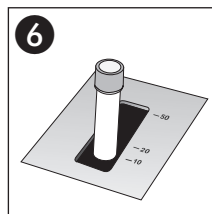

6
Messung
Measurement
Mesure
Medición