

# CSB 10000

## COD/DCO/DQO

# Test 0-23

REF 985 023



Rundküvettentest/ Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

1.00–10.00 g/L (1000–10000 mg/L)


620 nm

Method(e) / Método


0231

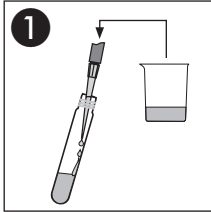


**1.0 mL** Probe  
Sample  
Echantillon  
Muestra

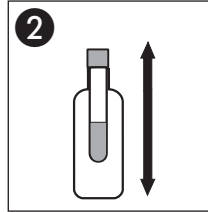


Sicherheitsgefäß  
Safety bottle  
Réceptient de sécurité  
Recipiente de seguridad

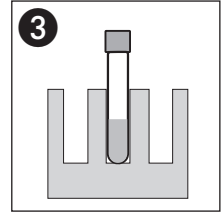


**148 °C / 2 h**

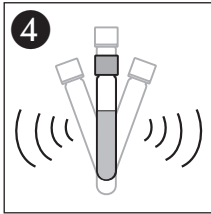


Schütteln  
Shake  
Agiter  
Agitar

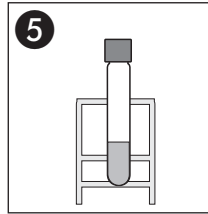

Abkühlen lassen  
Allow to cool  
Laisser refroidir  
Dejar enfriar

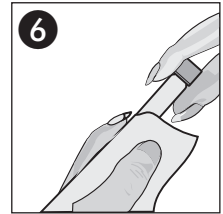

Säubern  
Clean  
Nettoyer  
Limpiar

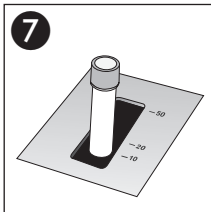

Messung  
Measurement  
Mesure  
Medición