

Chromat 5

Chromate/Cromato

Test 0-24

REF 985 024




Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

0.05–2.0 mg/L Cr(VI)


0.1–4.0 mg/L CrO₄²⁻

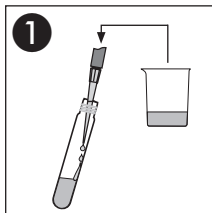


1
4.0 mL Probe
Sample
Echantillon
Muestra

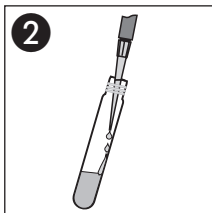


2
200 µL R2

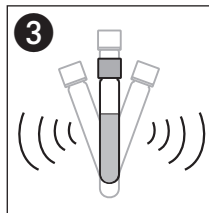


3
Schütteln
Shake
Agiter
Agitar

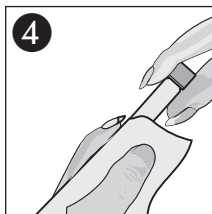

4
Säubern
Clean
Nettoyer
Limpiar

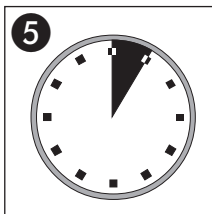

5
5'00 min

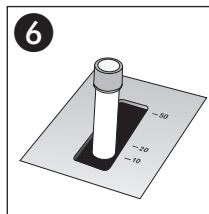


6
Messung
Measurement
Mesure
Medición

0.005–0.50 mg/L Cr(VI)

0.01–1.00 mg/L CrO₄²⁻

siehe Gebrauchsanweisung / see instructions for use / voir mode d'emploi / ver instrucciones de uso

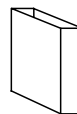

(50 mm)

Method(e) / Método

1241

1242