

# CSB 40

## COD/DCO/DQO

# Test 0-27

REF 985 027




Rundküvettentest/ Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica


2–40 mg/L

345 nm


Method(e) / Método  
0271

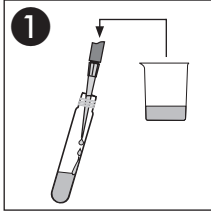


**2.0 mL** Probe  
Sample  
Echantillon  
Muestra

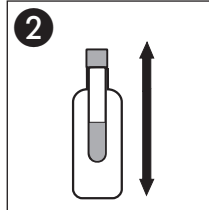


Sicherheitsgefäß  
Safety bottle  
Réceptient de sécurité  
Recipiente de seguridad

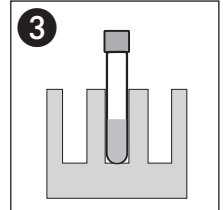


**148 °C / 2 h**

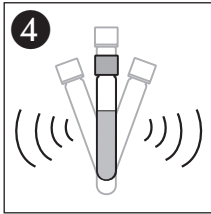


Schütteln  
Shake  
Agiter  
Agitar

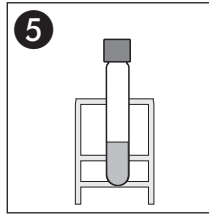

Abkühlen lassen  
Allow to cool  
Laisser refroidir  
Dejar enfriar

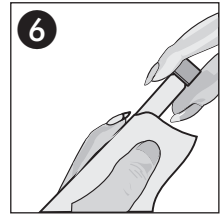

Säubern  
Clean  
Nettoyer  
Limpiar

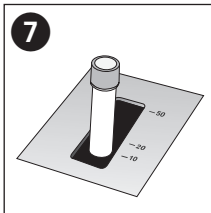

Messung  
Measurement  
Mesure  
Medición