

CSB 15000

COD/DCO/DQO

Test 0-28


REF 985 028




Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

1.0–15.0 g/L (1000–15000 mg/L)

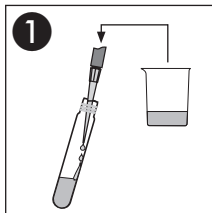

200 µL Probe
Sample
Echantillon
Muestra

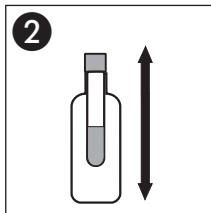


Sicherheitsgefäß
Safety bottle
Récipient de sécurité
Recipiente de seguridad


620 nm

Method(e) / Método


0281

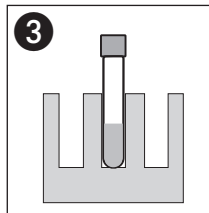


148 °C / 2 h

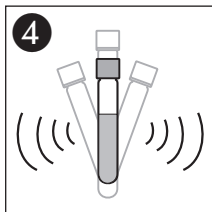


Schütteln
Shake
Agiter
Agitar

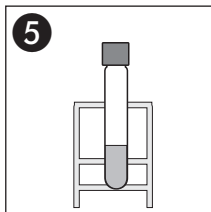

Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

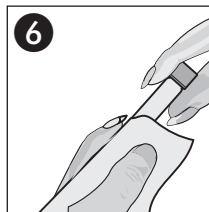

Säubern
Clean
Nettoyer
Limpiar

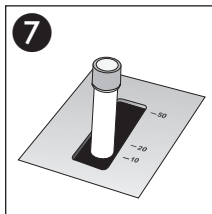

Messung
Measurement
Mesure
Medición