

# CSB 1500

## COD/DCO/DQO

# Test 0-29


REF 985 029




### Rundküvettentest/ Tube test

### Test en cuves rondes/ Prueba en cubeta cilíndrica

100–1500 mg/L

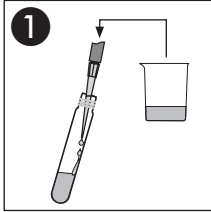

**2.0 mL** Probe  
Sample  
Echantillon  
Muestra

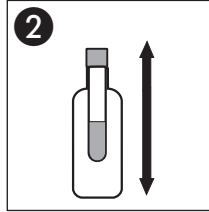


Sicherheitsgefäß  
Safety bottle  
Récipient de sécurité  
Recipiente de seguridad


620 nm

Method(e) / Método


**0291**

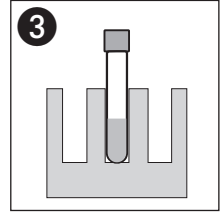


**148 °C/2 h**

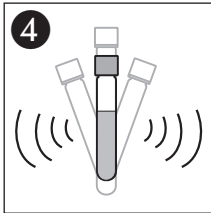


Schütteln  
Shake  
Agiter  
Agitar

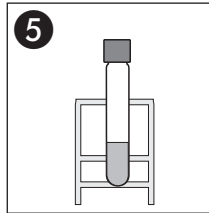

Abkühlen lassen  
Allow to cool  
Laisser refroidir  
Dejar enfriar

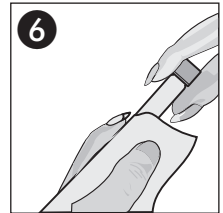

Säubern  
Clean  
Nettoyer  
Limpiar

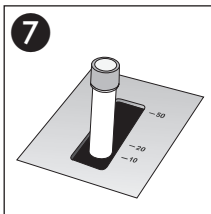

Messung  
Measurement  
Mesure  
Medición