

# CSB 600

## COD/DCO/DQO

### Test 0-30

REF 985 030



Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

50–600 mg/L


620 nm

Method(e) / Método


0301

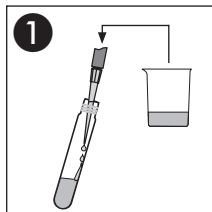


2.0 mL Probe  
Sample  
Echantillon  
Muestra

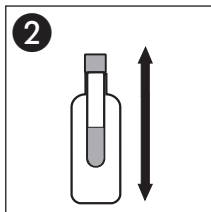


Sicherheitsgefäß  
Safety bottle  
Récipient de sécurité  
Recipiente de seguridad

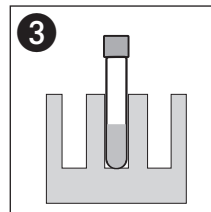


148 °C / 2 h

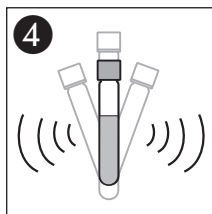


Schütteln  
Shake  
Agiter  
Agitar

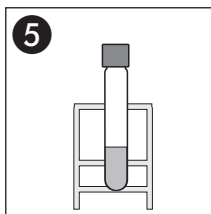

Abkühlen lassen  
Allow to cool  
Laisser refroidir  
Dejar enfriar

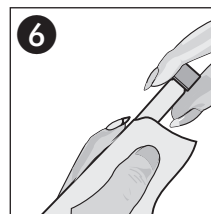

Säubern  
Clean  
Nettoyer  
Limpiar

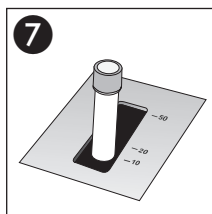

Messung  
Measurement  
Mesure  
Medición