

DEHA 1

Test 0-35

REF 985 035



Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

0.05–1.00 mg/L DEHA


540 nm

Method(e) / Método


0351

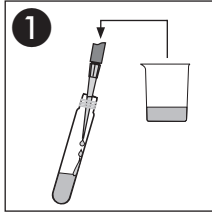


4.0 mL Probe
Sample
Echantillon
Muestra

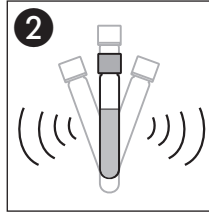


Schütteln
Shake
Agiter
Agitar

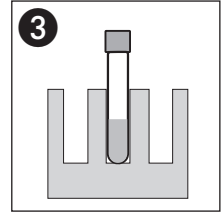


100 °C / 15 min

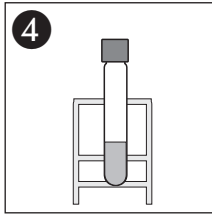


Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

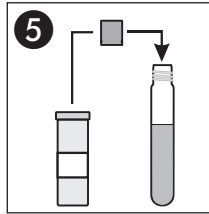


1 x NANOFIX R2

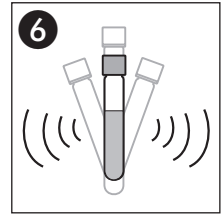


Schütteln
Shake
Agiter
Agitar

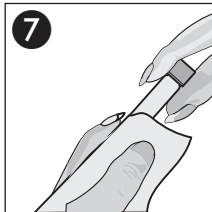

Säubern
Clean
Nettoyoy
Limpiar

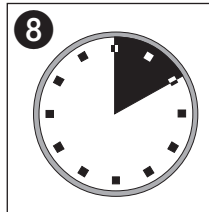

10'00 min

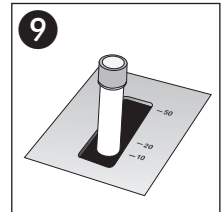

Messung
Measurement
Mesure
Medición