

# Fluorid 2

## Fluoride/Fluorure/Fluoruro

**Test 0-40**

REF 985 040



Rundküvettentest/ Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

0.1–2.0 mg/L F<sup>-</sup>


620 nm

Method(e) / Método


**0401**

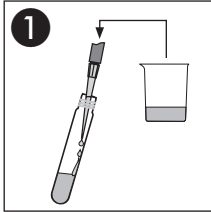


**2.0 mL** Probe  
Sample  
Echantillon  
Muestra

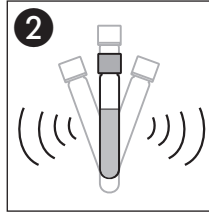


Schütteln  
Shake  
Agiter  
Agitar

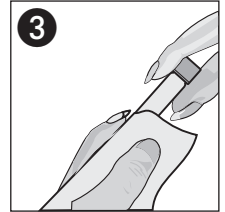

Säubern  
Clean  
Nettoyer  
Limpiar

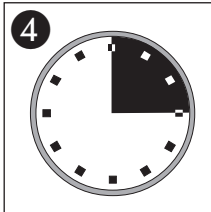

**15'00 min**

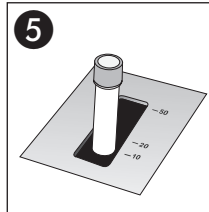

Messung  
Measurement  
Mesure  
Medición