

# Formaldehyd 8

## Formaldehyde/Formaldehydo

**Test 0-41**

REF 985 041



Rundküvettest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

0.1–8.0 mg/L HCHO


585 nm

Method(e) / Método


**0411**

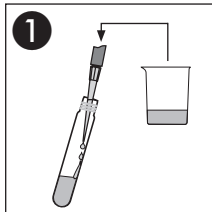


**2.0 mL Probe**  
Sample  
Echantillon  
Muestra

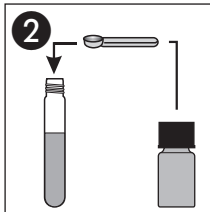


**1 x R2**

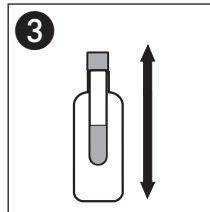


Sicherheitsgefäß  
Safety bottle  
Récipient de sécurité  
Recipiente de seguridad

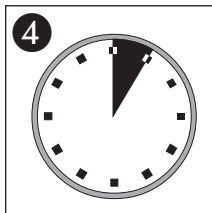

**5'00 min**

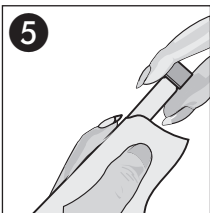

Säubern  
Clean  
Nettoyer  
Limpiar

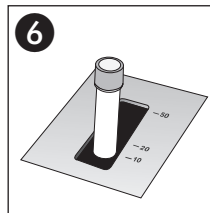

Messung  
Measurement  
Mesure  
Medición