

Härte 20

Hardness / Dureté / Durezza

Test 0-43

REF 985 043



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

1.0–20.0 °d / 0.2–3.6 mmol/L

1.0–25.0 °e / 2.0–36.0 °f


Gesamthärte / Total hardness / Dureté total / Durezza total:

540 nm


Method(e) / Método

0431/2


0433/4

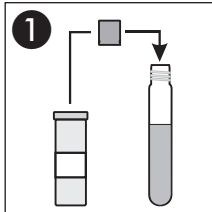


1 x **NANOFIX R2**

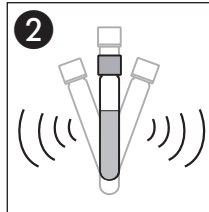


Kräftig schütteln
Shake well
Bien agiter
Agitar intensamente

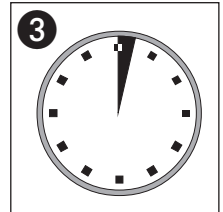


2'00 min

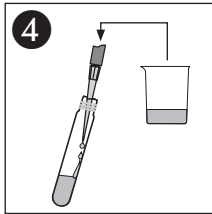


200 µL Probe
Sample
Echantillon
Muestra

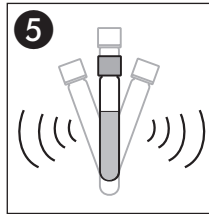


Schütteln
Shake
Agiter
Agitar

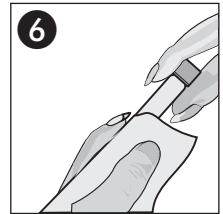


Säubern
Clean
Nettoyer
Limpiar

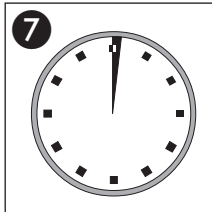

1'00 min

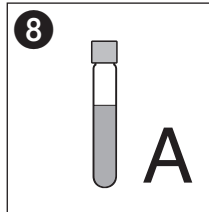

Lösung A
Solution A
Solution A
Solución A

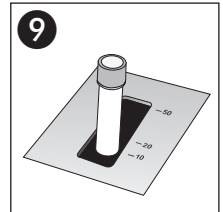

Messung
Measurement
Mesure
Medición