

Formaldehyd 10

Formaldehyde / Formaldehído

Test 0-46

REF 985 046



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

412 nm


Method(e) / Método

0.20–10.00 mg/L HCHO


0.20–10.00 ppm HCHO

0461


0462

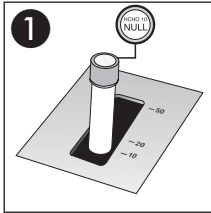


1
NULL messen
Measure blank
Lecture blanc
Lectura blanco

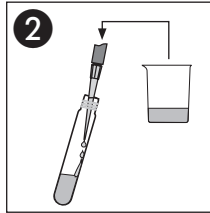


2
2.0 mL Probe
Sample
Echantillon
Muestra

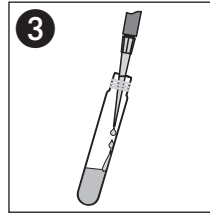


3
1.0 mL R2

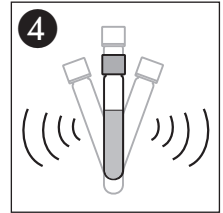


4
Schütteln
Shake
Agiter
Agitar

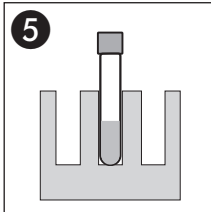


5
40 °C / 30 min

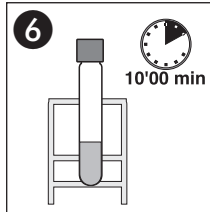


6
Abkühlen lassen
Allow to cool
Laisser refroidir
Dejar enfriar

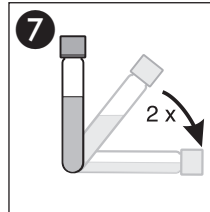


7
Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

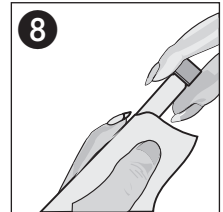

8
Säubern
Clean
Nettoyer
Limpiar

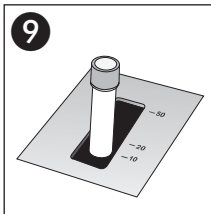

9
Messung
Measurement
Mesure
Medición

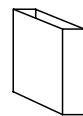

(50 mm)

Method(e) / Método

0.02–1.00 mg/L HCHO

1461

siehe Gebrauchsanweisung / see instructions for use /
voir mode d'emploi / ver instrucciones de uso