

Nichtionische Tenside 15

Nonionic surfactants / Tensio-actifs non ioniques / Tensioactivos no iónicos

Rundküvettestest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

0.3–15.0 mg/L Triton® X-100

Test 0-47


REF 985 047




620 nm

Method(e) / Método


0471

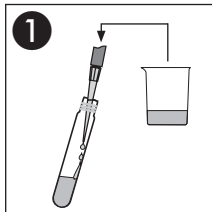


1
4.0 mL Probe
Sample
Echantillon
Muestra

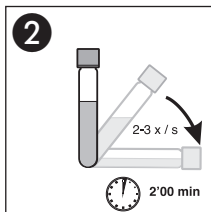


2
Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

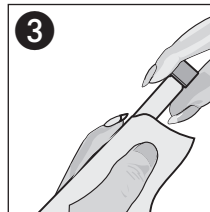

3
Säubern
Clean
Nettoyer
Limpiar

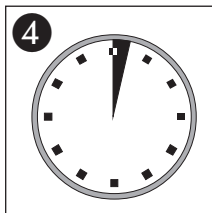

4
2'00 min
Phasentrennung abwarten
wait for phase separation
laisser séparer les phases
esperar la separación de fases

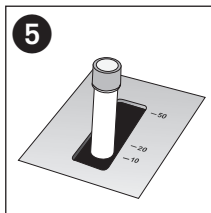

5
Messung
Measurement
Mesure
Medición