


Silber 3


Silver / Argent / Plata

Rundküvettest / Tube test
Test en cuves rondes / Prueba en cubeta cilíndrica


0.20–3.00 mg/L Ag⁺

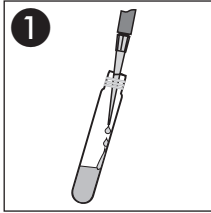


500 µL R2

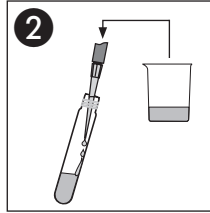


4.0 mL Probe
Sample
Echantillon
Muestra

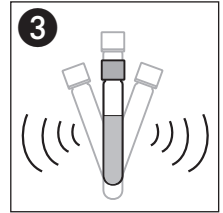


Schütteln
Shake
Agiter
Agitar

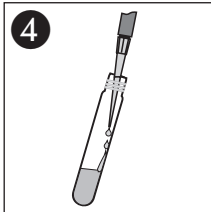


500 µL R3

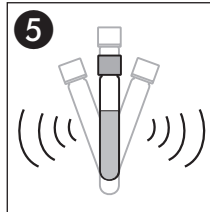


Schütteln
Shake
Agiter
Agitar

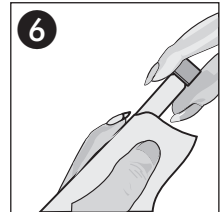


Säubern
Clean
Nettoyer
Limpiar

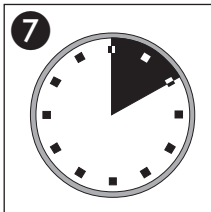

10'00 min

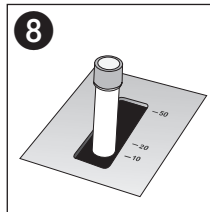

Messung
Measurement
Mesure
Medición

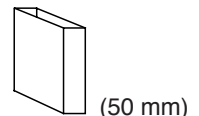

Method(e) / Método

0.08–0.50 mg/L Ag⁺

1491

siehe Sondervorschrift / see special instructions / voir procédures analytiques spéciales / ver instrucciones específicas de análisis