

# Org. Komplexbildner 10

## Complexing agents / Complexants organiques

### Agente formador de complejos org.

**Test 0-52**

REF 985 052



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica


540 nm

Method(e) / Método


0521

0.5–10.0 mg/L  $I_{\text{Bik}}$


Nullwert / Blanc value / Zéro / Cero

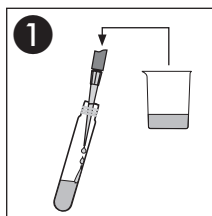


**1**  
4.0 mL dest. Wasser  
dist. water  
d'eau distillée  
agua dest.

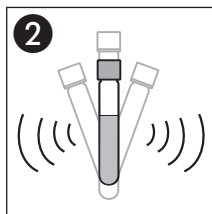


**2**  
Schütteln  
Shake  
Agiter  
Agitar

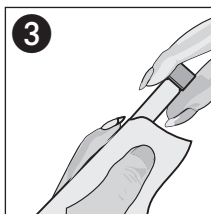

**3**  
Säubern  
Clean  
Nettoyer  
Limpiar

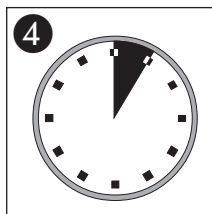


**4**  
5'00 min

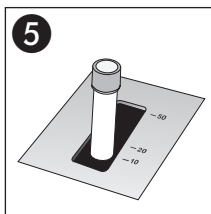


**5**  
Messung  
Measurement  
Mesure  
Medición


Messwert / Sample / Echantillon / Muestra

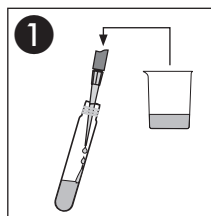


**1**  
4.0 mL Probe  
Sample  
Echantillon  
Muestra

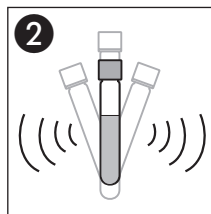


**2**  
Schütteln  
Shake  
Agiter  
Agitar

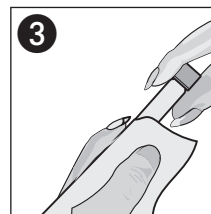

**3**  
Säubern  
Clean  
Nettoyer  
Limpiar

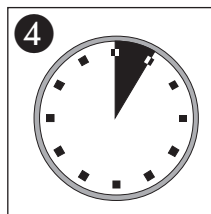

**4**  
5'00 min

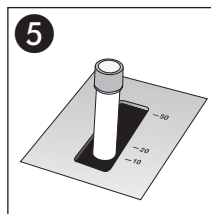

**5**  
Messung  
Measurement  
Mesure  
Medición