

# Kupfer 7

## Copper / Cuivre / Cobre

**Test 0-54**

REF 985 054



Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

585 nm


Method(e) / Método

**0541**


0.10–7.00 mg/L Cu<sup>2+</sup>

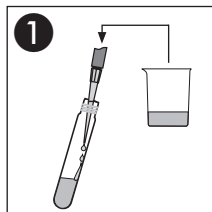


**1**  
4.0 mL Probe  
Sample  
Echantillon  
Muestra

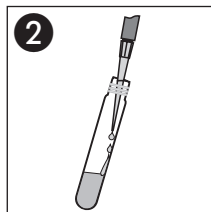


**2**  
200 µL R2

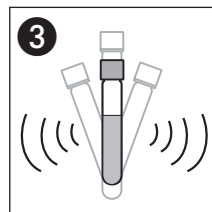


**3**  
Schütteln  
Shake  
Agiter  
Agitar

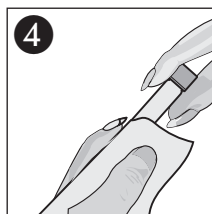

**4**  
Säubern  
Clean  
Nettoyer  
Limpiar

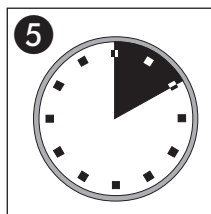

**5**  
10'00 min

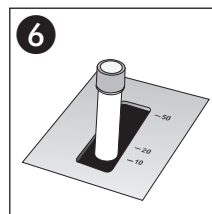

**6**  
Messung  
Measurement  
Mesure  
Medición