

ortho-Phosphat 45

ortho Phosphate / orto Fosfato

Test 0-55

REF 985 055



Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

5.0–50.0 mg/L PO₄-P

15–150 mg/L PO₄³⁻

10–115 mg/L P₂O₅


690 nm

Method(e) / Método


0554

0555


0556

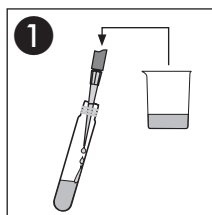


200 µL Probe
Sample
Echantillon
Muestra

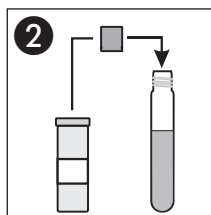


1 x NANOFIX R3

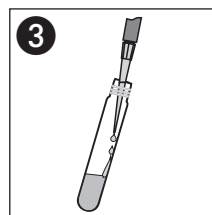


200 µL R4

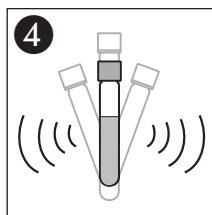


Schütteln
Shake
Agiter
Agitar

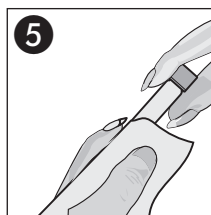

Säubern
Clean
Nettoyer
Limpiar

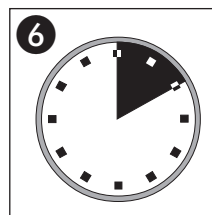

10'00 min

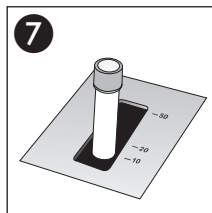

Messung
Measurement
Mesure
Medición