

Mangan 10

Manganese / Manganèse / Manganeso

Test 0-58

REF 985 058



Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

470 nm


Method(e) / Método

0.1–10.0 mg/L Mn


0581

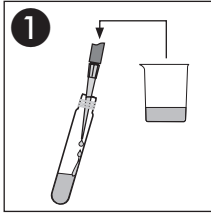


1
4.0 mL Probe
Sample
Echantillon
Muestra

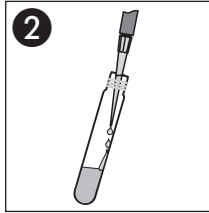


2
500 μ L R2

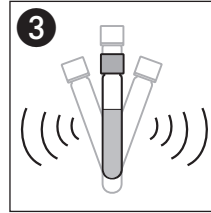


3
Schütteln
Shake
Agiter
Agitar

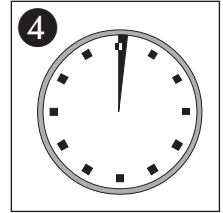


4
1'00 min

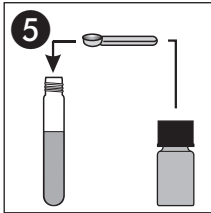


5
1 x R3

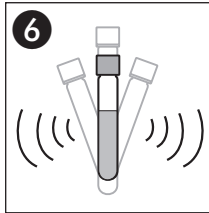


6
Kräftig schütteln
Shake well
Bien agiter
Agitar intensamente

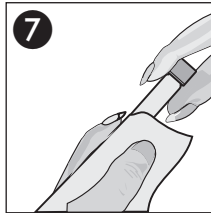


7
Säubern
Clean
Nettoyer
Limpiar

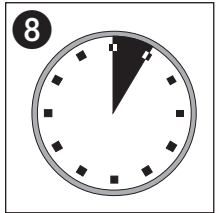

8
5'00 min

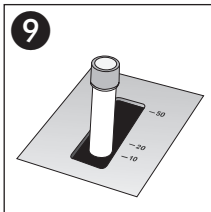

9
Messung
Measurement
Measure
Medición

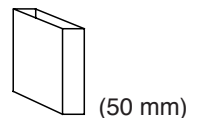

Method(e) / Método

0.02–2.00 mg/L Mn

1581

siehe Gebrauchsanweisung / see instructions for use /
voir mode d'emploi / ver instrucciones de uso