

# Nitrat 8

## Nitrate / Nitrato

**Test 0-65**

REF 985 065




Rundküvettest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

0.30–8.00 mg/L  $\text{NO}_3\text{-N}$


1.3–35.0 mg/L  $\text{NO}_3^-$

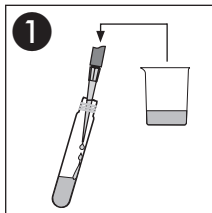


**1**  
500  $\mu\text{L}$  Probe  
Sample  
Echantillon  
Muestra

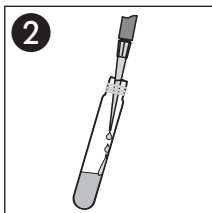


**2**  
500  $\mu\text{L}$  R2

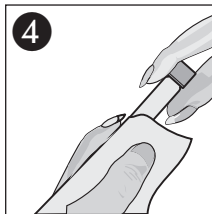


**4**  
Säubern  
Clean  
Nettoyer  
Limpiar

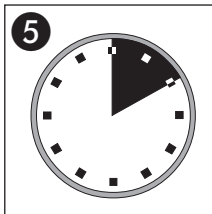

**5**  
10'00 min

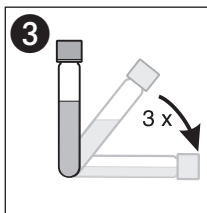

**3**  
Umschwenken  
Shake gently  
Secouer légèrement  
Mezclar volteándolo

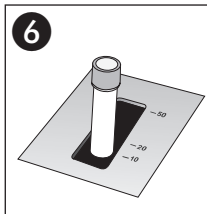

**6**  
Messung  
Measurement  
Mesure  
Medición