

Nitrat 250

Nitrate / Nitrato

Test 0-66


REF 985 066




Rundküvettestest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica


4–60 mg/L $\text{NO}_3\text{-N}$
20–250 mg/L NO_3^-

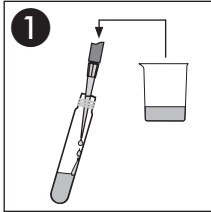


200 μL Probe
Sample
Echantillon
Muestra

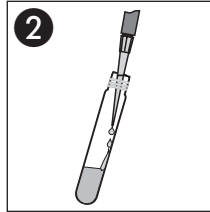


500 μL R2

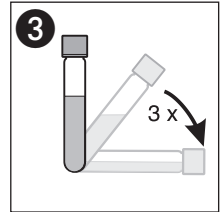


Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

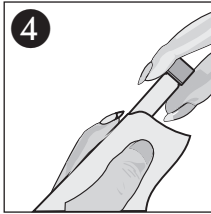

Säubern
Clean
Nettoyer
Limpiar

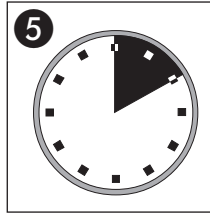

10'00 min

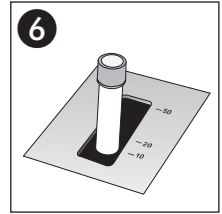

Messung
Measurement
Mesure
Medición