

Nitrit 2

Nitrite / Nitrito

Test 0-68

REF 985 068



Rundküvettest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

0.003–0.460 mg/L NO₂-N


0.02–1.50 mg/L NO₂⁻

540 nm


Method(e) / Método

0681


0682

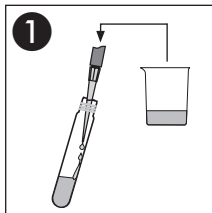


1
4.0 mL Probe
Sample
Echantillon
Muestra

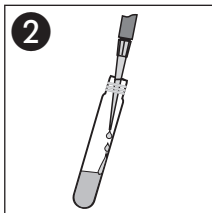


2
200 µL R2

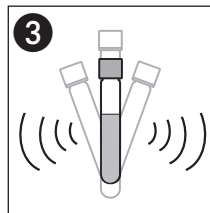


3
Schütteln
Shake
Agiter
Agitar

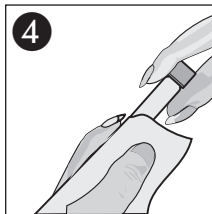

4
Säubern
Clean
Nettoyer
Limpiar

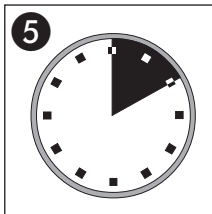

5
10'00 min

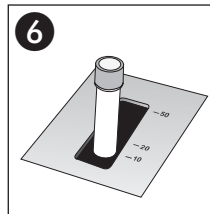

6
Messung
Measurement
Mesure
Medición