

Nitrit 4

Nitrite / Nitrito

Test 0-69

REF 985 069



Rundküvettestest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

0.1–4.0 mg/L NO₂-N


0.3–13.0 mg/L NO₂⁻

540 nm


Method(e) / Método

0691


0692

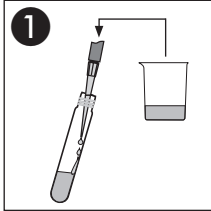


200 µL Probe
Sample
Echantillon
Muestra

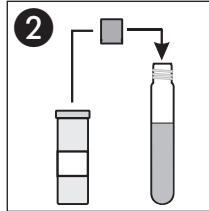


1 x NANOFIX R2

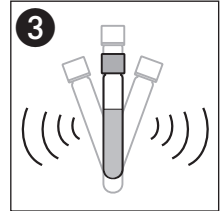


Schütteln
Shake
Agiter
Agitar

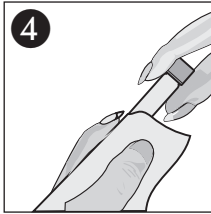

Säubern
Clean
Nettoyer
Limpiar

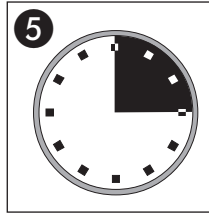

15'00 min

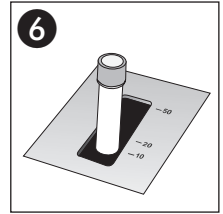

Messung
Measurement
Mesure
Medición