

POC 200

Polyoxycarbonsäuren / Polyoxycarboxylic acids

Acides polyoxycarboxyliques / Acidos polioxicarboxílicos

Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

Test 0-70

REF 985 070



436 nm


Method(e) / Método

0704


2–40 mg/L Polystabil[®] KWI

Anwendung im Kühlwasser / Application for cooling water


Application dans l'eau de circulation / Aplicación en agua de circulación

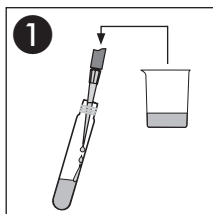


1
4.0 mL Probe
Sample
Echantillon
Muestra

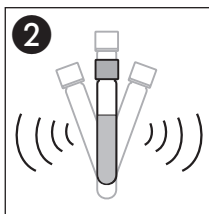


2
Schütteln
Shake
Agiter
Agitar

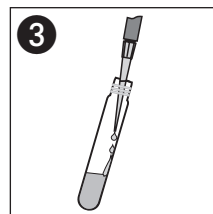


3
1.0 mL R2

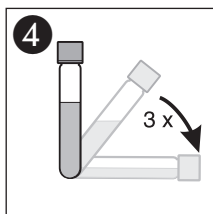


4
Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

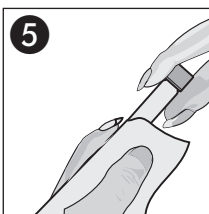


5
Säubern
Clean
Nettoyer
Limpiar

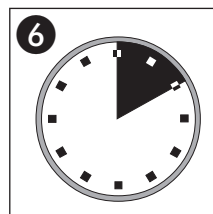

6
10'00 min

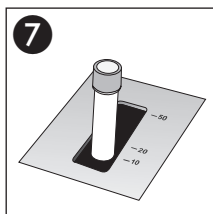

7
Messung
Measurement
Mesure
Medición

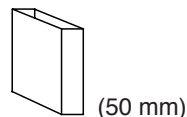

(50 mm)

Method(e) / Método

1.0–10.0 mg/L Polystabil[®] KWI

1704

siehe Gebrauchsanweisung / see instructions for use /
voir mode d'emploi / ver instrucciones de uso