

Sulfid 3

Sulphide/Sulfure/Sulfuro

Test 0-73

REF 985 073



Rundküvettentest/ Tube test


Test en cuves rondes/ Prueba en cubeta cilíndrica

0.05–3.00 mg/L S²⁻


620 nm

Method(e) / Método


0731

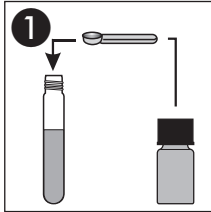


1 x R2

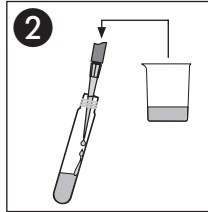


4.0 mL Probe
Sample
Echantillon
Muestra

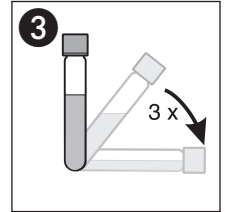


Umschwenken
Shake gently
Secouer légèrement
Mezclar volteándolo

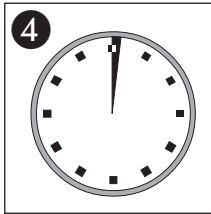


1'00 min

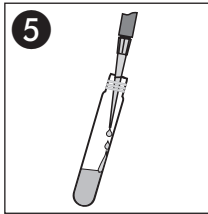


200 µL R3

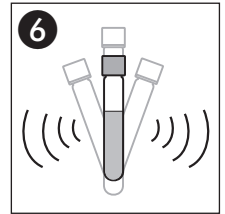


Schütteln
Shake
Agiter
Agitar

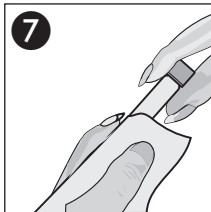

Säubern
Clean
Nettoyer
Limpiar

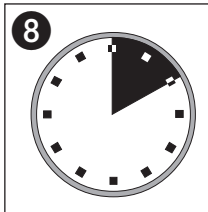

10'00 min

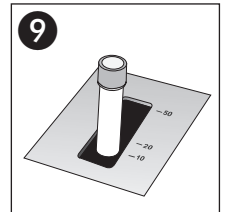

Messung
Measurement
Mesure
Medición