


Phenol-Index 5

Phenolic index / Indice phénol / Indice fenolico


Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica


0.2–5.0 mg/L

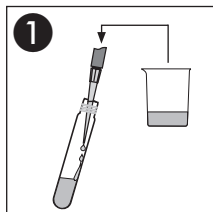


4.0 mL Probe
Sample
Echantillon
Muestra

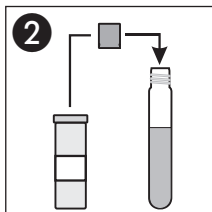


1 x **NANOFIX R2**

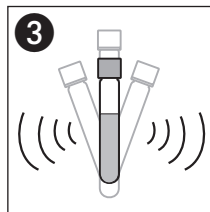


Schütteln
Shake
Agiter
Agitar

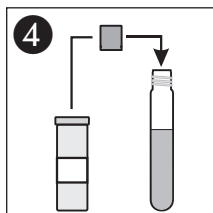


1 x **NANOFIX R3**

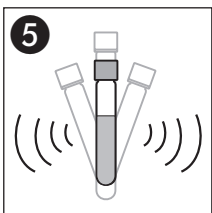


Schütteln
Shake
Agiter
Agitar

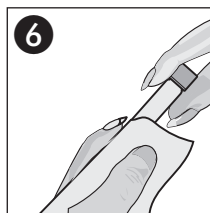

Säubern
Clean
Nettoyer
Limpiar

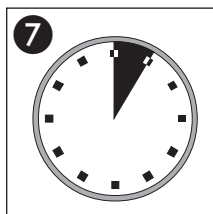

5'00 min

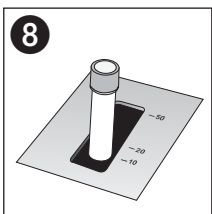

Messung
Measurement
Mesure
Medición

Test 0-74

REF 985 074



520 nm

Method(e) / Método

0741