

Phenol-Index 5

Phenolic index / Indice phénol / Indice fenolico

Test 0-74

REF 985 074
+ 918 929



Rundküvettentest / Tube test
Test en cuves rondes / Prueba en cubeta cilíndrica


470 nm

Method(e) / Método


1742

0.2–5.0 mg/L


Nullwert / Blanc value / Zéro / Cero

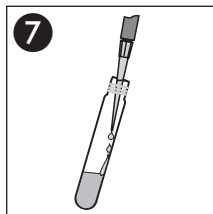


4.0 mL MIBK

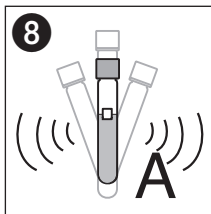


30 s kräftig schütteln
shake well
bien agiter
agitar intensamente

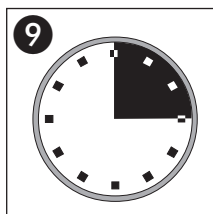

15'00 min

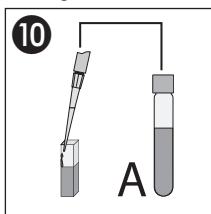


3 mL in Küvette
into cell
dans cuve
en cubeta

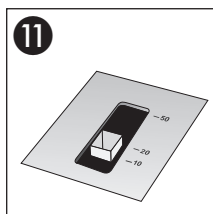


Messung
Measurement
Mesure
Medición


Messwert / Sample / Echantillon / Muestra

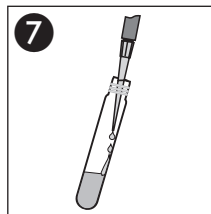


4.0 mL MIBK

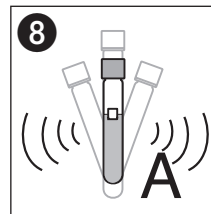


30 s kräftig schütteln
shake well
bien agiter
agitar intensamente

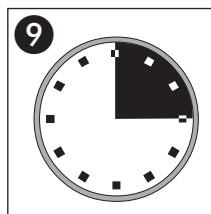

15'00 min

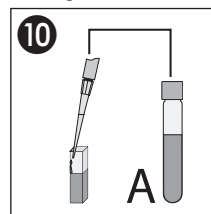

3 mL in Küvette
into cell
dans cuve
en cubeta

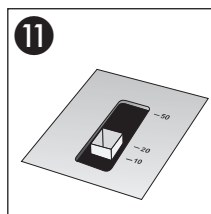

Messung
Measurement
Mesure
Medición