

ortho-Phosphat 50

ortho Phosphate / orto Fosfato

Test 0-79

REF 985 079



Rundküvettest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

10.0–50.0 mg/L $\text{PO}_4\text{-P}$

30–150 mg/L PO_4^{3-}

23–110 mg/L P_2O_5


436 nm

Method(e) / Método


0794

0795


0796

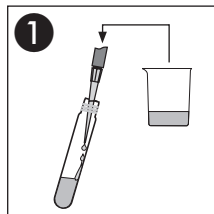


1
4.0 mL Probe
Sample
Echantillon
Muestra

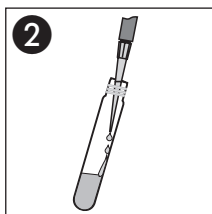


2
1 mL R3

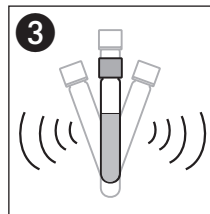

3
Schütteln
Shake
Agiter
Agitar

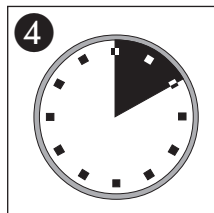

4
10'00 min

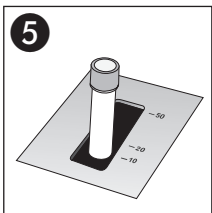

5
Messung
Measurement
Mesure
Medición