

Resthärte 1

Residual Hardness / Dureté résiduelle

Dureza residual

Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica

Test 0-84

REF 985 084



540 nm


Method(e) / Método

0841 / 2


0843 / 4

0.02–1.00 °d / 0.004–0.180 mmol/L


0.03–1.25 °e / 0.04–1.78 °f

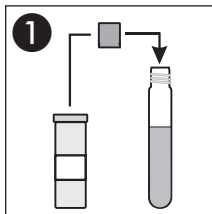


1 x **NANOFIX R2**

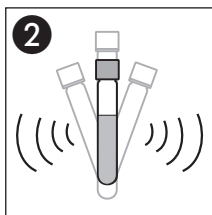


Kräftig schütteln
Shake well
Bien agiter
Agitar intensamente

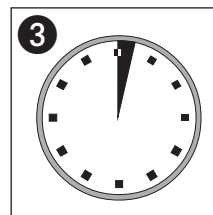


2'00 min

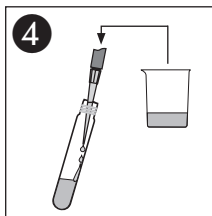


5.0 mL Probe
Sample
Echantillon
Muestra

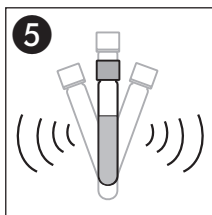


Schütteln
Shake
Agiter
Agitar

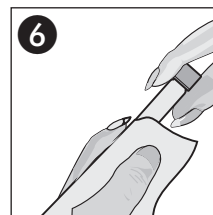

Säubern
Clean
Nettoyer
Limpiar

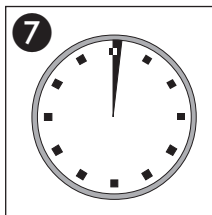

1'00 min

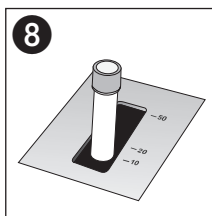

Messung
Measurement
Mesure
Medición