

Sulfat 200

Sulphate/Sulfate/Sulfato

Test 0-86

REF 985 086



Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

10–200 mg/L SO_4^{2-}


436 nm

Method(e) / Método


0861

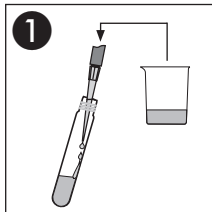


1
4.0 mL Probe
Sample
Echantillon
Muestra

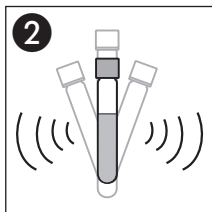


2
Schütteln
Shake
Agiter
Agitar

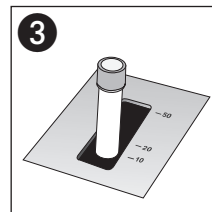


3
Nullwert
Blank value
Zéro
Cero

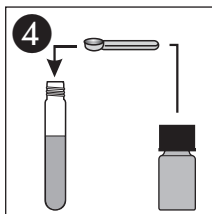


4
1 x R2

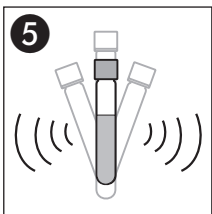


5
10 s Schütteln
Shake
Agiter
Agitar

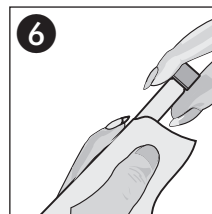

6
Säubern
Clean
Nettoyer
Limpiar

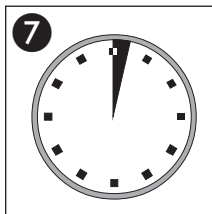

7
2'00 min

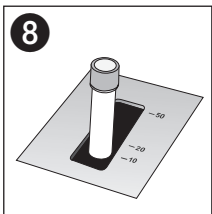

8
Messung
Measurement
Mesure
Medición