

# Sulfat 1000

## Sulphate / Sulfate / Sulfato

**Test 0-87**

REF 985 087



Rundküvettentest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

436 nm


Method(e) / Método

0871


200–1000 mg/L  $\text{SO}_4^{2-}$

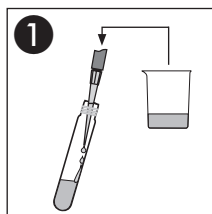


**1.0 mL** Probe  
Sample  
Echantillon  
Muestra

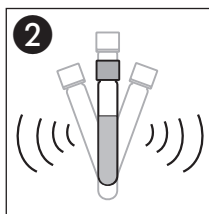


**Schütteln**  
Shake  
Agiter  
Agitar

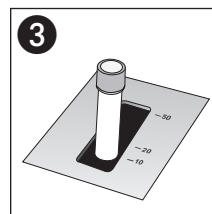


**Nullwert**  
Blank value  
Zéro  
Cero

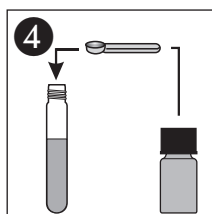


**1 x R2**

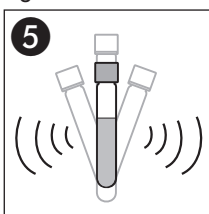


**10 s Schütteln**  
Shake  
Agiter  
Agitar

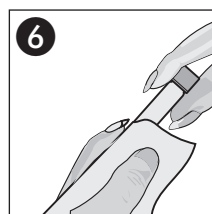

**Säubern**  
Clean  
Nettoyer  
Limpiar

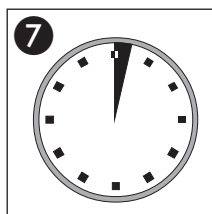

**2'00 min**

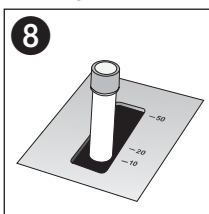

**Messung**  
Measurement  
Mesure  
Medición