

Sulfit 10

Sulphite/Sulfite/Sulfito

Test 0-89


REF 985 089




Rundküvettentest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica


0.2–10.0 mg/L SO_3^{2-}

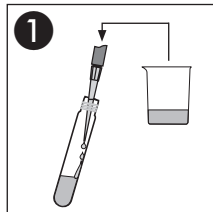


1
4.0 mL Probe
Sample
Echantillon
Muestra

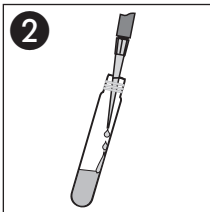


2
200 μL R2

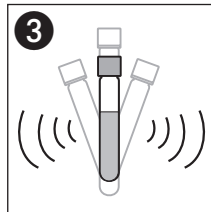


3
Schütteln
Shake
Agiter
Agitar

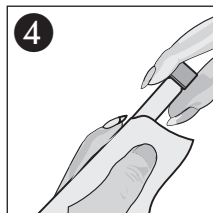

4
Säubern
Clean
Nettoyer
Limpiar

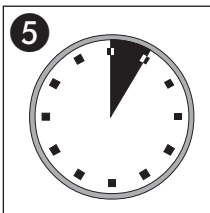


5
5'00 min

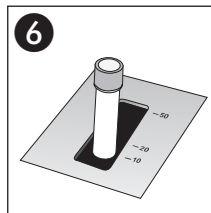

6
Messung
Measurement
Mesure
Medición

0.05–2.40 mg/L SO_3^{2-}

siehe Gebrauchsanweisung / see instructions for use / voir mode d'emploi / ver instrucciones de uso

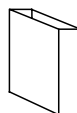

(50 mm)

Method(e) / Método

1891