

# Sulfit 100

Sulphite/Sulfite/Sulfito

**Test 0-90**

REF 985 090



Rundküvettest/ Tube test  
Test en cuves rondes/ Prueba en cubeta cilíndrica

470 nm

Method(e) / Método


0901

5–100 mg/L  $\text{SO}_3^{2-}$


Nullwert/ Blanc value/ Zéro/ Cero




**Leere** Rundküvette  
**Empty** test tube  
Cuve ronde **vide**  
Tubo de test **vacío**

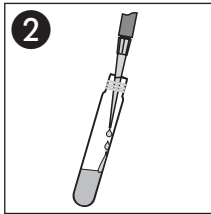

1.0 mL R1

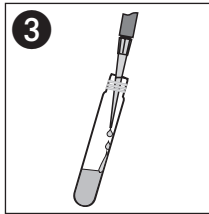


200  $\mu\text{L}$  R2

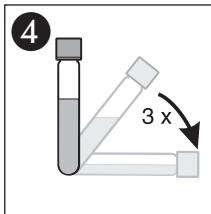


Umschwenken  
Shake gently  
Secouer légèrement  
Mezclar volteándolo


Messwert/Sample/Echantillon/Muestra

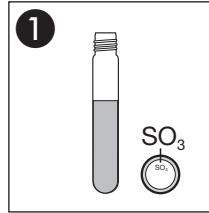

**Sulfit**-Rundküvette  
Test tube **Sulphite**  
Cuve ronde **Sulfite**  
Tubo de test **Sulfito**

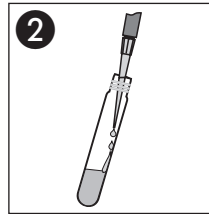

200  $\mu\text{L}$  R2

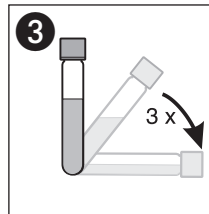

Umschwenken  
Shake gently  
Secouer légèrement  
Mezclar volteándolo