

Thiocyanat 50

Thiocyanate / Tiocianato

Test 0-91

REF 985 091



Rundküvettestest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

0.5–50.0 mg/L SCN⁻


470 nm

Method(e) / Método


0911

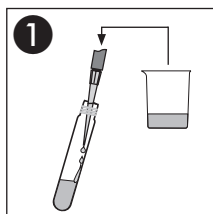


1
4.0 mL Probe
Sample
Echantillon
Muestra

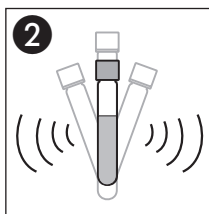

2
Schütteln
Shake
Agiter
Agitar

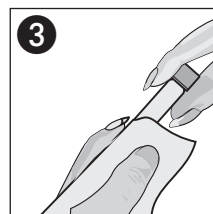

3
Säubern
Clean
Nettoyer
Limpiar

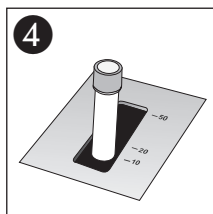

4
Messung
Measurement
Measure
Medición