

Zink 4

Zinc

Test 0-96

REF 985 096



Rundküvettest / Tube test


Test en cuves rondes / Prueba en cubeta cilíndrica

620 nm


Method(e) / Método

0961


0.10–4.00 mg/L Zn²⁺

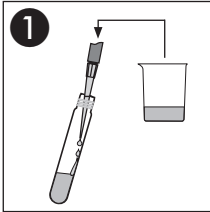


1
4.0 mL Probe
Sample
Echantillon
Muestra

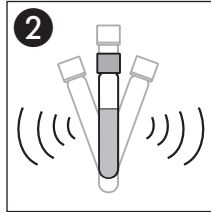


2
Schütteln
Shake
Agiter
Agitar

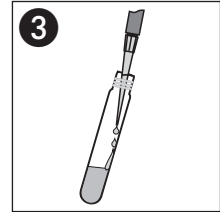


3
200 µL R2

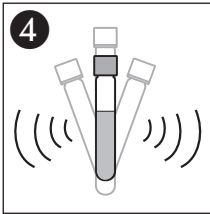


4
Schütteln
Shake
Agiter
Agitar

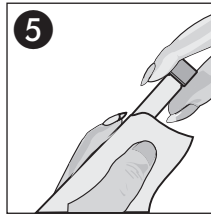

5
Säubern
Clean
Nettoyer
Limpiar

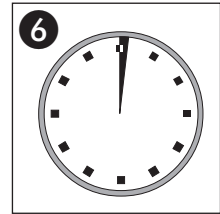

6
1'00 min

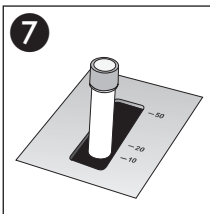

7
Mesung
Measurement
Mesure
Medición