

Aluminium 07

Aluminio

Test 0-98

REF 985 098



Rundküvettentest / Tube test


540 nm

Test en cuves rondes / Prueba en cubeta cilíndrica


Method(e) / Método

0.02–0.70 mg/L Al³⁺



0981

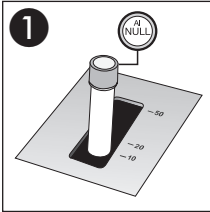


1
NULL messen
Measure blank
Lecture blanc
Lectura blanco

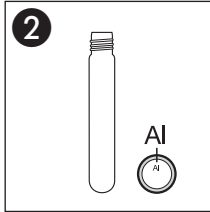


2
Aluminium-Rundküvette 500 µL R2
Test tube **Aluminium**
Cuve ronde **Aluminium**
Tubo de test **Aluminio**

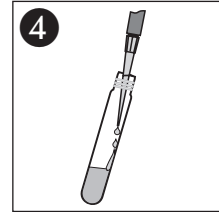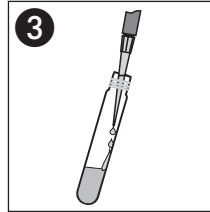


3
500 µL R3

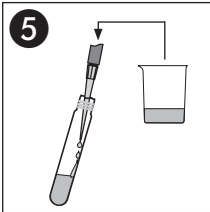


5
4.0 mL Probe
Sample
Echantillon
Muestra

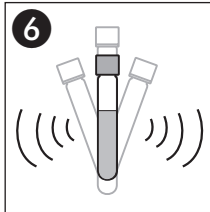


6
Schütteln
Shake
Agiter
Agitar

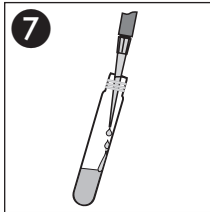


7
500 µL R4

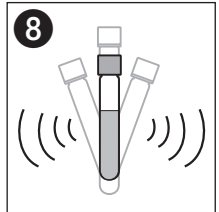


8
Schütteln
Shake
Agiter
Agitar

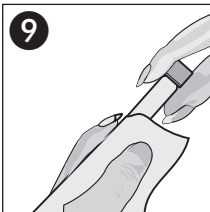

9
Säubern
Clean
Nettoyer
Limpiar

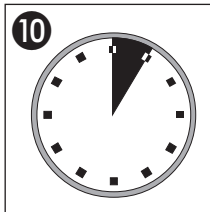

10
5'00 min

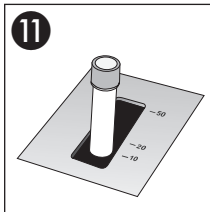

11
Messung
Measurement
Mesure
Medición