

TOC 600 COT

Test 0-99


REF 985 099




Rundküvettest / Tube test

Test en cuves rondes / Prueba en cubeta cilíndrica


40–600 mg/L C

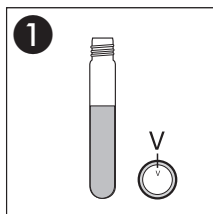


Verdünnungsküvette
Dilution test tube
Cuve ronde de dilution
Tubo de dilución

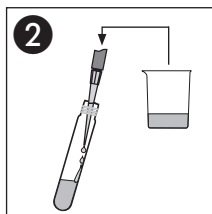


1.0 mL Probe
Sample
Echantillon
Muestra

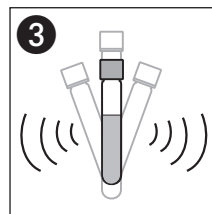


Schütteln
Shake
Agiter
Agitar

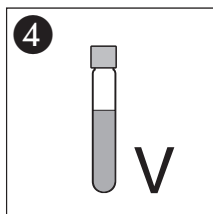


Verdünnungslösung V
Diluted solution V
Solution de dilution V
Solución de dilución V

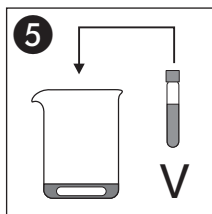


Verdünnungslösung
Diluted solution
Solution de dilution
Solución de dilución

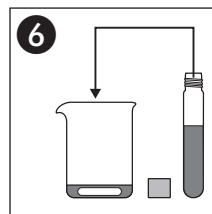


500 µL R0

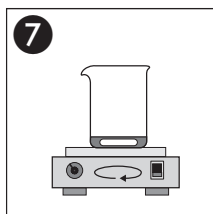

10 min Rühren
Stirr
Agiter
Agitar

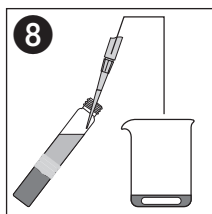

4.0 mL in TOC-Küvette
into tube TOC
dans cuve COT
en tubo TOC

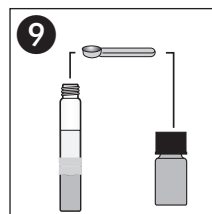

1 x R2

585 nm

Method(e) / Método

0991