

Ethanol 1000

Etanol

Test 8-38

REF 985 838




Rundküvettestest / Tube test
Test en cuves rondes / Prueba en cubeta cilíndrica

620 nm


Method(e) / Método

0.1–1.0 g/L EtOH
0.013–0.13 Vol% EtOH


8381
8382

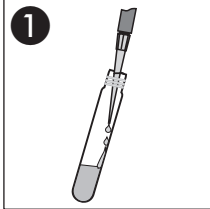


4.0 mL R1

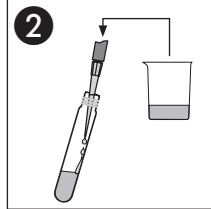


0.5 mL Probe
Sample
Echantillon
Muestra

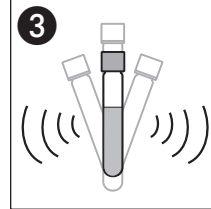


Schütteln
Shake
Agitar
Agitar

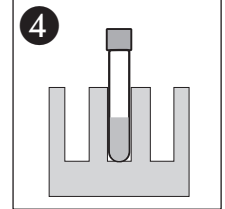


25 °C / 20 min

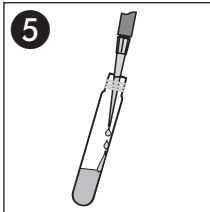


100 µL R2

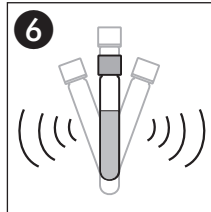


Schütteln
Shake
Agitar
Agitar

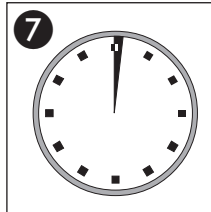


1'00 min

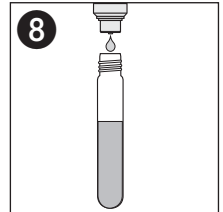


2 x R3

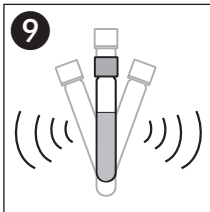


Schütteln
Shake
Agitar
Agitar

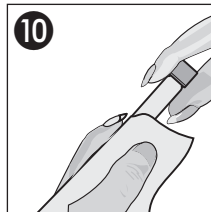

Säubern
Clean
Nettoyer
Limpiar

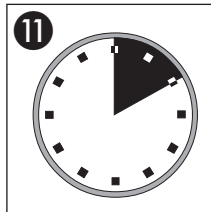

10'00 min

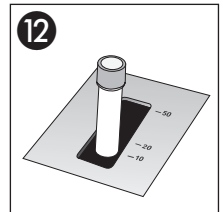

Messung
Measurement
Mesure
Medición