

Methanol 15

Méthanol / Metanol

Test 8-59

REF 985 859




Rundküvettest / Tube test
Test en cuves rondes / Prueba en cubeta cilíndrica

620 nm


Method(e) / Método

8591


0.2–15.0 mg/L MeOH

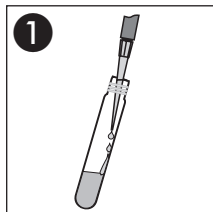


3.0 mL R1

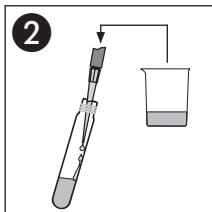


1.5 mL Probe
Sample
Echantillon
Muestra

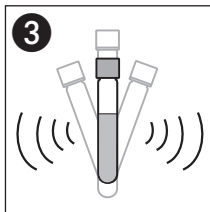


Schütteln
Shake
Agiter
Agitar

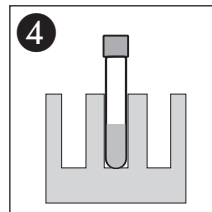


25 °C / 30 min

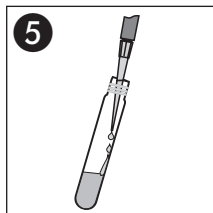


100 µL R2

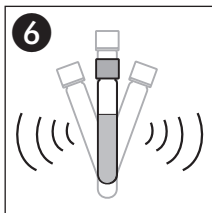


Schütteln
Shake
Agiter
Agitar

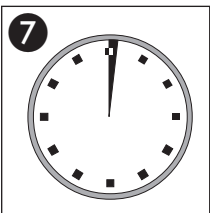


1'00 min

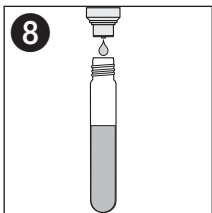


2 x R3

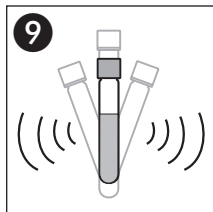


Schütteln
Shake
Agiter
Agitar

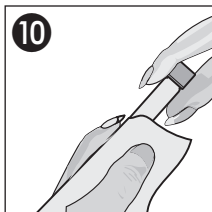

Säubern
Clean
Nettoyer
Limpiar

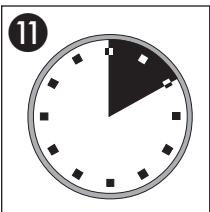

10'00 min

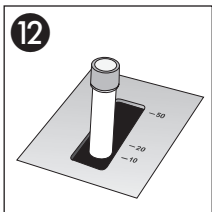

Messung
Measurement
Mesure
Medición