


visicolor® ECO Aluminium

Aluminio

Messbereich / Range / Domaine de mesure / Rango
0.10 - 0.50 mg/l Al³⁺


REF 931 006
931 206

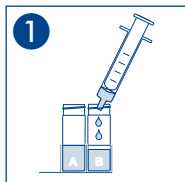


1
2 x 5 ml Probe
Sample
Echantillon
Muestra




2
1 Al-1

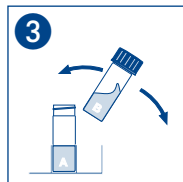


3
Schütteln
Shake
Agiter
Agitar

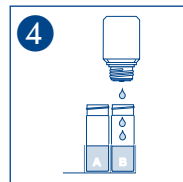


4 Al-2

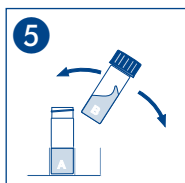


5
Schütteln
Shake
Agiter
Agitar

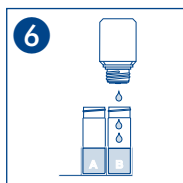


6 Al-3

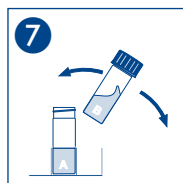


7
Schütteln
Shake
Agiter
Agitar

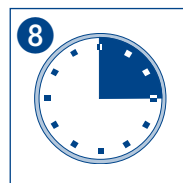

8
15'00 min

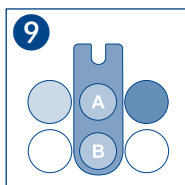

9
Messung
Measurement
Mesure
Medición

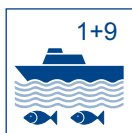

1+9
Meerwasser
Sea water
Eau de mer
Agua de mar