


visicolor® ECO Ammonium 3 Amonio




Messbereich / Range / Domaine de mesure / Rango
0.2 - 3 mg/l NH_4^+


REF 931 008
931 208

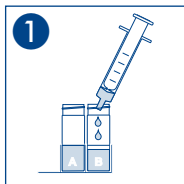

1
2 x 5 ml Probe
Sample
Echantillon
Muestra

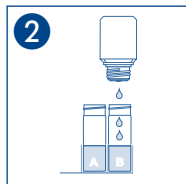


2
10 NH_4 -1

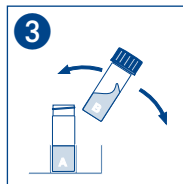


3
Schütteln
Shake
Agiter
Agitar




4
1 NH_4 -2

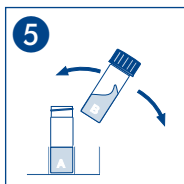


5
Schütteln
Shake
Agiter
Agitar

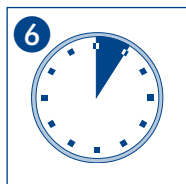


6
5'00 min

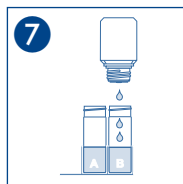


7
4 NH_4 -3

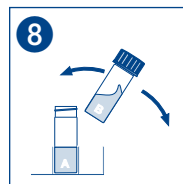


8
Schütteln
Shake
Agiter
Agitar

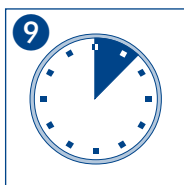

9
7'00 min

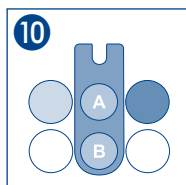

10
Messung
Measurement
Mesure
Medición

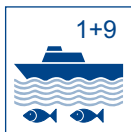

1+9

Meerwasser
Sea water
Eau de mer
Agua de mar